



Proximel[®] and Proximel[®] Ag

Non-Border Non-Border

Expect more from your silicone foam dressing.

Proximel[®] Non-Border and **Proximel[®] Non-Border Ag** are silicone foam dressings that feature a unique tapered edge. This unique design allows the dressing to lay smooth around wound edges, helping prevent exudate leakage.¹ This makes it an ideal choice for management of moderate to high exuding wounds, as well as for use underneath compression.² The silicone adhesive on Proximel[®] Non-Border and Proximel[®] Non-Border Ag helps the dressing stay in place during bandage wrapping, is gentle on the skin and is easily removed.

Proximel[®] Non-Border Ag offers antimicrobial activity against a wide range of pathogens³ and provides a sustained silver release over time, allowing the dressing to be worn for up to 7 days.

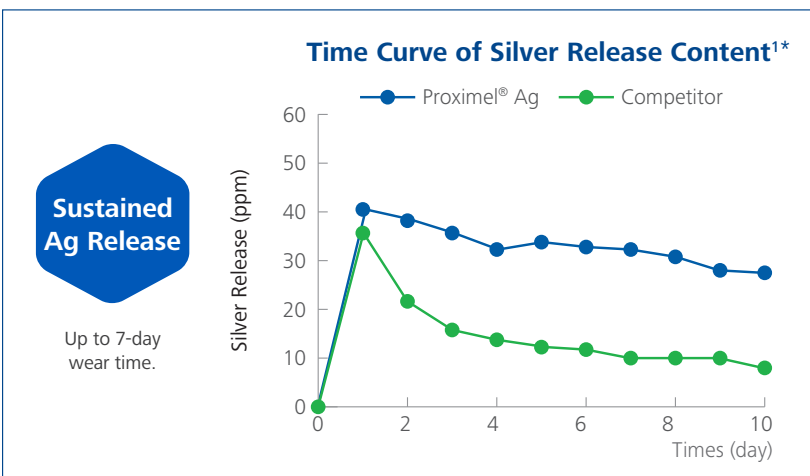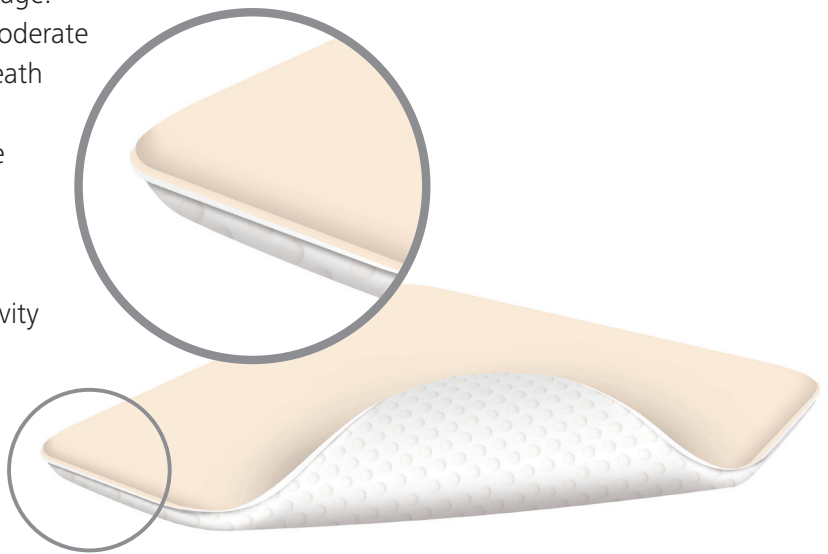


Features and Benefits:

- Pain-free removal and easy to reposition
- Tapered edge design is ideal for use under compression
- Helps reduce the risk of maceration
- Fast exudate absorption
- Dynamic fluid management provides optimal moist wound environment

*Applies to Proximel[®] Non-Border Ag only.

Proximel® Non-Border: Ideal for use under compression²

Tapered edge locks in moisture¹ while laying gently around wounds

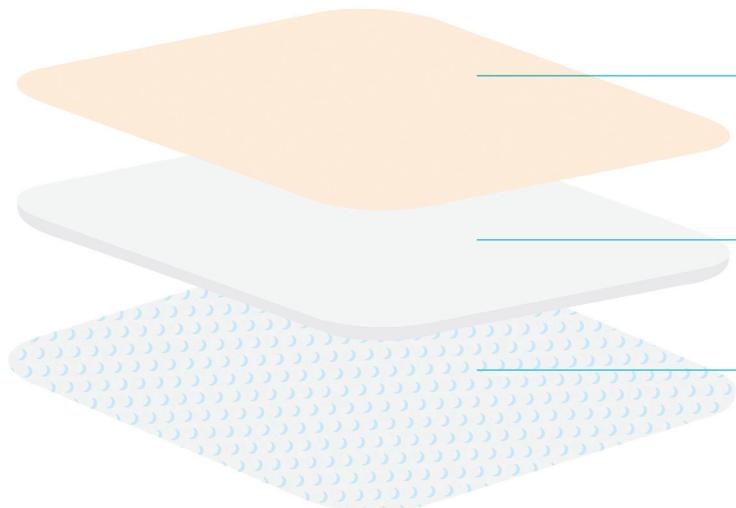

1 Protective Backing

High MVTR
Smooth and comfortable
Waterproof and bacteria barrier

2 Absorbent Polyurethane Foam

Great fluid handling capability
Offers additional comfort and cushioning protection

3 Silicone Wound Contact Layer

Unique silicone oil and perforation technology
Distributes wound exudate fast and evenly
Gentle and secure with pain free removal
Does not strip away epidermal cells



Sustained silver release ensures effectiveness against key pathogens.*



Proximel's distinctive properties allow you to confidently **lift** and **reposition** the dressing while reducing the risk of adherence failure.



High MVTR rate supports good moisture balance.



Superabsorbent core supports **fast absorption** and helps minimize maceration risk.

Proximel® Non-Border

Order number	Dressing size	Dressings per box
15713100	4" x 4" / 10cm x 10cm	0
15713101	4" x 5" / 10cm x 12cm	10
15713102	4" x 8" / 10cm x 20cm	5
15713103	6" x 6" / 15cm x 15cm	5
15713104	8" x 8" / 20cm x 20cm	5
15713105	Heel, 5" x 8" / 13cm x 20cm	5

Proximel® Non-Border Ag

Order number	Description	Dressings per box
15713117	4" x 4" / 10cm x 10cm	10
15713118	6" x 6" / 15cm x 15cm	5
15713119	8" x 8" / 20cm x 20cm	5

*Applies to Proximel® Non-Border Ag only. ¹ Data on file. ² Data on file. ³ Data on file.

Contact your local HARTMANN representative, or visit us at www.hartmannusa.com for more information.

HARTMANN USA, Inc.
Rock Hill, SC 29730
1-800-243-2294

www.hartmannusa.com

XLIT 2794 Rev. 4 (0719)



Going further
for health